



NE1BD-14



Designed for applications that require continuous boiling. This bath also has the flexibility where the digital control can be regulated from a gentle simmer to a vigorous boil.

Comprising of a stainless steel tank housed in a durable paint finished outer case. Painted surfaces feature a tough anti-bacterial finish with good chemical resistance. Complete with a removable stainless steel perforated false base, drain tap and clear lift off lid.


Fitted with a sheathed immersion element. With the heater directly immersed in water it provides rapid boiling. A standard feature is a constant level device, which maintains a fixed level of liquid thereby eliminating the need to regularly top up the bath. Not fitted on NE1BD-18W.

NE1BD - 14

NE1BD - 18W

NE1BD - 28

GENERAL			
Chamber Capacity	14 Litres	18 Litres	28 Litres
Temperature Range	Simmer 80°C to boiling at 105°C		
Temperature Control	Digital, LED display showing actual or set, resolution 1°C		
Heater Power (220-240V)	2500W		
Voltage	230V / 50-60Hz		
Water Level Alarm	Features low liquid level float device, heating automatically suspended		
Sterilising Tray	<input checked="" type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>
Clear Lift Off Lid Supplied		<input checked="" type="checkbox"/>	
Drain Tap Fitted		<input checked="" type="checkbox"/>	
Constant Level Device	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
DIMENSIONS			
External (mm)	362w x 332d x 290h		538w x 332d x 290h
Internal (mm)	325w x 300d x 150h	325w x 300d x 200h	500w x 300d x 200h

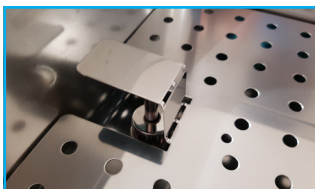

Low liquid float device

A low liquid float device protects the heaters and a red warning light on the front panel alerts users to refill the bath.



NE1BD-18W

Raise and lower sterilising tray along with pair of guide rail assemblies comes as standard. Provides a safe method of lifting sterilised items from boiling water. Excess hot water drains away.

NE1BD Digital Boiling Baths

Accessories - Lids




Stainless Steel Gable Lid

Bath	NE1BD-14	NE1BD-18W	NE1BD-28
Accessory Code	N/A		SL1-22SS



Stainless Steel Flat Lid

Bath	NE1BD-14	NE1BD-18W	NE1BD-28
Accessory Code	LD-14	N/A	LD-22

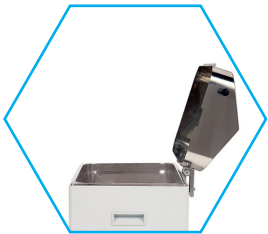

Stainless Steel Hinged Gable Lid

Bath	NE1BD-14	NE1BD-18W	NE1BD-28
Accessory Code	SL1-14H/SS		N/A



Clear Hinged Gable Lid

Bath	NE1BD-14	NE1BD-18W	NE1BD-28
Accessory Code	GL1-14		GL1-22



Stainless Steel Concentric Ringed Lid

Bath	NE1BD-14	NE1BD-18W	NE1BD-28
Accessory Code	SLR1-14	N/A	SLR1-22
	SLR2-14		

SLR1-14: 4 x 114mm holes

SLR2-14: 6 x 83mm holes

SLR1-22: 6 x 114mm holes



Polypropylene Spheres (Pack of 200)

Bath	NE1BD-14	NE1BD-18W	NE1BD-28
No. of packs	2 Packs		4 Packs
Accessory Code	BP0368		

NE1BD Digital Boiling Baths

Accessories



Stainless Steel Test Tube Racks

Bath	NE1BD-14	NE1BD-18W	NE1BD-28
No. of Racks	3		5
Accessory Code	6870 : 26 holes x 17mm DIA		
	6871 : 16 holes x 26mm DIA		
	6872 : 36 holes x 13mm DIA		
	6873 : 18 holes x 19mm DIA		
	6900 : 12 holes x 32mm DIA		

Universal Accessories

CAL	Calibration Record
-----	--------------------